

# samaritan™ PAD

## Public Access Defibrillator



HeartSine®

### Compact, easy-to-use, lifesaving technology.

Sudden Cardiac Arrest (SCA) is a leading cause of death globally. Response time is critical for survival. HeartSine's samaritan PAD was designed especially for use in public areas by providing a sophisticated defibrillator for adult or pediatric use, inside a lightweight and easy-to-operate system.




### Compact in size, long on ability.

**Compact.** The samaritan PAD is lighter (2.4 lbs.) and smaller than other defibrillators.

**Durable.** The samaritan PAD resists shock and vibration and carries an IP56 Rating, the industry's highest rating against dust and water. It can be taken and used virtually anywhere, even in the most inclement conditions. It also carries a 7-year unit warranty.

**Advanced technology.** SCOPE™\* Biphasic technology automatically optimizes energy output for each person, and has been reported to provide significantly better performance in removing ventricular fibrillation (VF) by the third shock.\*\*

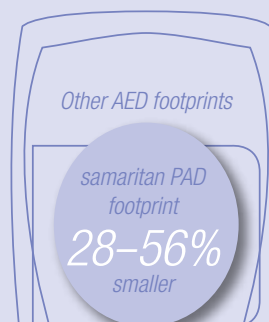

2.4 lbs. light.

### Easy-to-follow visual and verbal guides.

**User-friendly.** The samaritan PAD features easy-to-understand visual and oral prompts that guide a user through the process.

**Two-button operation.** Only two buttons, ON and SHOCK, are required, providing straightforward operation.

**Always ready.** A System Status Ready Indicator flashes to show that the complete system is operational and ready for use. Device automatically runs self check each week.



Visual cues prompt pad placement



Stand clear of the patient



It is safe to touch the patient


### Real economy for the real world.

**Two parts, one expiration date.** Pad-Pak™ cartridge combines battery and electrode pads, with one expiration date to monitor.

**Low cost of ownership.** Cartridge typically has a shelf-life of 3.5 years from date of manufacture, offering significant savings over other defibrillators that require separate battery and pad units.



Pad-Pak™ and Pediatric-Pak™, with pre-attached electrodes

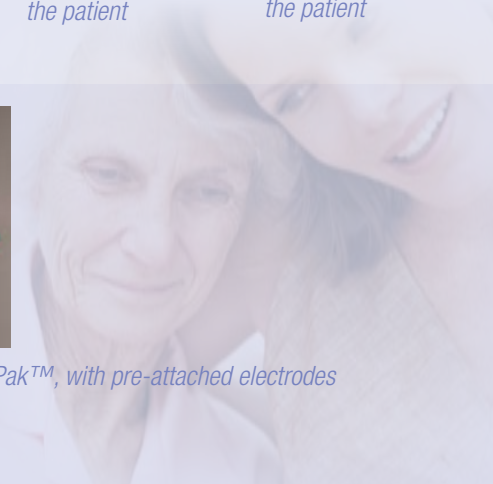



# Technical Overview

## Physical

**Size:** 8.0 in x 7.25 in x 1.9 in (20cm x 18.4cm x 4.8cm)  
**Weight:** 2.4 lbs (1,1 kg) including Pad-Pak™ Battery

## Defibrillator

**Waveform:** SCOPE® (self-compensating output pulse envelope) Biphasic waveform. Optimized biphasic waveform compensates energy, slope and envelope for patient impedance

## Energy Selection

Pre-configured factory settings for escalating energy are AHA/ERC 2005  
**Adult**  
1. Shock 150J 2. Shock 150J 3. Shock 200J  
**Paediatric**  
1. Shock 50J 2. Shock 50J 3. Shock 50J

## Patient Analysis System

**Method:** Evaluates patient's ECG, signal quality, electrode contact integrity and patient impedance to determine if defibrillation is required  
**Sensitivity/Specificity:** Meets AAMI DF80:2003 requirements

## Environmental

**Operating Temperature:** 0°C to 50°C (+32°F to +122°F)  
**Storage Temperature:** -10°C to 50°C (+14°F to +122°F)  
**Relative Humidity:** 5% to 95% (non-condensing)  
**Water Resistance:** IEC 60529/EN 60529 IP56  
**Altitude:** 0 to 15,000 feet (0 – 4,575 meters)  
**Shock:** MIL STD 810F Method 516.5, Procedure I (40G's)  
**Vibration:** MIL STD 810F Method 514.5+  
Category 4 Truck Transportation – US Highways  
Category 7 Aircraft – Jet 737 & General Aviation (Exposure)  
**EMC:** EN 60601-1-2, Second Edition: 2001  
**Radiated Emissions:** CISPR11: 1997 +1A:1999+A2:2001 Group 1 Class B  
**Electrostatic Discharge Immunity:** EN61000-4-2:1995 (8KV)+A1:1998+A2:2000  
**RF Immunity:** EN61000-4-3:1996, 80 MHz – 2.5 GHz, (10V/m)+A1:1998+A2:2000  
**Magnetic Field Immunity:** EN61000-4-8:1993 (3 A/m)+A1:2000  
**Aircraft:** RTCA / DO – 160D: 1997, Section 21 (Category M)  
**Falling height:** 1 meter

## Event Documentation

**Type:** Internal memory  
**Memory Capacity:** 45 minutes of ECG (full disclosure) and event/incident recording  
**Playback Capabilities:** Custom USB cable directly connected to PC and Saver™ EVO Windows-based data review software

## Materials Used

**PAD SAM300P:** ABS, Santoprene. Printed circuit board with electronic components.  
**PAD Cartridge:** Battery: Lithium Manganese Dioxide  
Housing: ABS – Electrodes: Hydrogel, Silver, Aluminium and Polyester

## Combined Electrode and Battery Cartridge

**Pad-Pak™**  
**Shelf Life:** Check Expiration Date (typically 3.5 years from manufacture)  
**Weight:** 0.44 lbs (0.2kg)  
**Size:** 3.93 in x 5.24 in x .94 in (10cm x 13.3cm x 2.4cm)  
**Battery Type:** Lithium Manganese Dioxide (LiMnO<sub>2</sub>)  
18V, 0.8 Amp Hrs  
>30 shocks at 200J or 6 hours of continuous monitoring.  
**Capacity:**

## Adult Electrodes (Pad-Pak-01)

**Standby Life:** Typically 3.5 years from manufacture date  
**Electrodes:** Samaritan® disposable defibrillation pads are supplied as standard with each device  
Anterior-lateral  
**Placement:** 100cm<sup>2</sup>  
**Active Gel Area:** 3.5 ft (1m)  
**Cable Length:**

## Children Electrodes (Pad-Pak-02)

**Standby Life:** Typically 3.5 years from manufacture date  
**Electrodes:** Samaritan® disposable defibrillation pads are supplied as standard with each device  
Anterior-Posterior  
**Placement:** 100cm<sup>2</sup>  
**Active Gel Area:** 3.5 ft (1m)cable  
**Cable Length:**



\*Self Compensating Output Pulse Envelope technology automatically compensates energy, slope and pulse envelope for the patient.

\*\*Efficacy of Distinct Energy Delivery Protocols Comparing Two Biphasic Defibrillators for Cardiac Arrest, Walsh, McClelland, Owens, Anderson, Turner, Adgey; *The American Journal of Cardiology*, Vol. 94, Aug.1, 2004

## Corporate Headquarters:

HeartSine Technologies, Inc.  
105 Terry Drive,  
Newtown, PA. 18940  
Tel: (215) 860 8100  
Fax: (215) 860 8192  
Web: heartsine.com

## Manufacturing Facilities

HeartSine Technologies, Ltd  
Canberra House, 203 Airport Road West  
Belfast, Northern Ireland BT3 9ED  
Tel: +44 (0) 28 9093 9400  
Fax: +44 (0) 28 9093 9401

## Distributed by:

AED Everywhere, Inc.  
3241 Nite Court  
Fort Collins, CO 80525  
877-751-5300



The products described in this brochure all meet the applicable European Medical Directive requirements.